IPS **Empress**[®] Esthetic Speed

Instructions for Use

Page 2

Gebrauchsinformation

Seite 5

Mode d'emploi

Page 8

Istruzioni d'uso

Pagina 11

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Página 14

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Página 17

Bruksanvisning

Sidan 20

Brugsanvisning

Side 23

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Sivu 26

Bruksanvisning

Side 29

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Pagina 32

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www.ivoclarvivadent.com

ivoclar
vivadent[®]
technical

The phosphate-bonded IPS Empress® Esthetic Speed investment material (Type 1, Class 2) is used for the pressing of IPS Empress Esthetic ingots in the Ivoclar Vivadent Programat press furnaces.

Indication

Investing models (wax, modelling resin) in preparation of the pressing procedure with IPS Empress Esthetic ingots.

Contraindication

All other applications not listed in the indications.

Delivery form

- IPS Empress Esthetic Speed Powder, 25 bags, 100 g each, including 1 measuring cup
- IPS Empress Esthetic Speed Liquid, 0.5 litre

Storage

- Storage at 12–28 °C / 54–82 °F
- Store powder in a dry place.
- The liquid must not be exposed to temperatures below +5°C / 41 °F (frost susceptible!). **If frozen once, the liquid is unusable (crystallization) and must not be used anymore.**

Processing temperatures

- Room temperature 18 °C – max. 23 °C / 64 °F – max. 73 °F
- Any other processing temperature decisively affects the setting behaviour.
- A higher processing temperature accelerates the setting behaviour.
- A lower processing temperature delays the setting behaviour.

Investment ring system and press plunger

The IPS Empress Esthetic investment ring system (ring base, ring gauge, silicone ring) and the Alox plunger/ One-Way-Plunger are coordinated with IPS Empress Esthetic. Components of different systems are not compatible. Use therefore only components of this system.

Mixing ratio 100 g powder : 27 ml diluted liquid

Indication	Concentration approx.	100-g investment ring Liquid : dist. water	200-g investment ring Liquid : dist. water
Inlay MO + OD	50%	13.5 ml : 13.5 ml	27 ml : 27 ml
Inlays MOD	60%	16 ml : 11 ml	32 ml : 22 ml
Veneers	50%	13.5 ml : 13.5 ml	27 ml : 27 ml
Crowns	60%	16 ml : 11 ml	32 ml : 22 ml

The indications regarding the liquid concentration are approximates and may be adjusted, if required. The total quantity of liquid (liquid + dist. water) must not be altered. Do not dilute the liquid to below 50%.

Processing

Excess of separating agents may cause reactions with the investment material and affect the press result. Disperse excess material with compressed air before investment.

- Fill liquid into the mixing cup.
- Add powder.
- Thoroughly mix investment material with a spatula for 20 seconds until even wetting of the powder is achieved.

Vacuum mixer

Vacuum mixing time (at room temperature)	1.5 min
Speed	350 rpm
Processing time	5.5 min
Setting time	min. 30 min – max. 45 min

The processing time depends on the material temperature, mixing quantity, mixing time and mixing intensity. A higher material temperature and/or a longer mixing time shortens the processing time.

Investment

- Mix max. 400 g investment material in one go.
- Closely observe mixing ratio!
- Investment is carried out on a shaker under slight vibration. Allow the investment ring to set without manipulating it.

Some aspects influencing the setting expansion of the investment material

- Temperature of the investment material and the liquid
- Residual water in the mixing cup, atmospheric humidity

Preheating / holding times

Starting temperature (Switch on the preheating furnace in time)	850 °C / 1562 °F
Holding time after reaching the final temperature of 850 °C / 1562 °F	100-g investment ring – min. 45 min
	200-g investment ring – min. 60 min

If more than 2 investment rings are preheated at the same time, the holding time at the final temperature has to be prolonged by 15 minutes per additional investment ring.

Note

If several Speed investments are to be conducted, they should be invested consecutively and placed into the preheating furnace at an interval of approx. 20 minutes. When placing the investment rings in the preheating furnace, make sure that the furnace temperature does not drop substantially. The stipulated holding time counts from the point when the preheating temperature has been reached again.

IPS Empress Esthetic	
Ingot	Preheating
Alox plunger	Preheating
One-Way-Plunger	No preheating

The IPS Unitrax ensures that the ingots and the Alox plungers are optimally preheated and positioned in the preheating furnaces.

Important information

- Do not use a debubbler (tenside). Reactions with the investment material may occur.
- Consistent results can only be achieved under consistent conditions and if the Instructions for Use are observed.
- Dilute liquid exclusively with distilled or de-ionized water.
- Liquid concentration: The more distilled or de-ionized water is added to the mixing liquid, the lower the setting expansion. Additionally, the strength of the investment material is reduced. Do not dilute the liquid to below 50%.
- Use only dry, clean instruments to process the investment material. Do not use plaster mixing cups.
- Do not invest under pressure since this changes the product properties.
- The setting time must be closely observed, since it decisively influences the expansion.
- Always remove pressed investment rings from the press furnace immediately after the end of the program and place them on a grid for cooling.

To ensure thorough heating of the investment rings in the preheating furnace, the following points have to be observed:

- Load the preheating furnace to maximally half of the available utility space.
- Always place the investment rings in the rear part of the preheating furnace.
- Place the investment rings as quickly as possible in the preheating furnace. Make sure that the furnace temperature does not drop substantially.
- Always place the investment rings in the preheating furnace with the opening facing down.
- If more than two 200-g investment rings have to be preheated in the preheating furnace, we recommend investing them consecutively and placing them in the preheating furnace at an interval of approximately 20 minutes.
- Do not place the investment rings in direct contact with each other. Allow air circulation.


IMPORTANT! Safety notes

- Do not open the furnace during the first 20 minutes of preheating: Burn hazard caused by wax fumes that might ignite in the air!
- This material contains quartz, which may cause lung diseases (silicosis or cancer) when inhaled. Prevent the inhalation of dust and wear a suitable dust mask.

Physical data

(according to ISO 15912:2006)

- | | | |
|-----------------------------|--------------|--------------------------|
| - Flow characteristics: | 110 mm | Liquid concentration 50% |
| - Initial set: | 10.5 minutes | Liquid concentration 50% |
| - Compressive strength: | 5 MPa | Liquid concentration 50% |
| - Linear thermal expansion: | 0.9 % | Liquid concentration 50% |

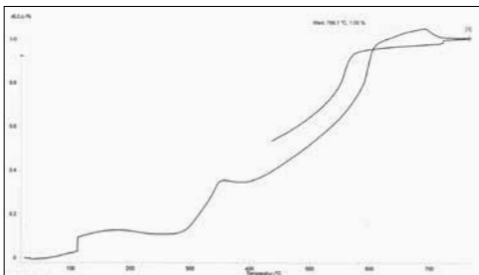This material has been developed solely for use in dentistry. Processing should be carried out strictly according to the Instructions for Use. Liability cannot be accepted for damages resulting from failure to observe the Instructions or the stipulated area of application. The user is responsible for testing the products for their suitability and use for any purpose not explicitly stated in the Instructions. Descriptions and data constitute no warranty of attributes and are not binding.

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Ivoclar Vivadent Limited
Ivoclar Vivadent, Inc.