

SR PHONARES®

IMPLANT ESTHETICS



Anterior Artistry

SR PHONARES® anterior teeth embody the essence of natural dentition. Each mould has been sculpted to deliver the ultimate lifelike result.


Outstanding Esthetics

- Unparalleled surface texture
- Natural translucency, opalescence and fluorescence



Functional Design

- Set and fit design ensures proximal harmony
- Characterized setups can be achieved for a distinct look

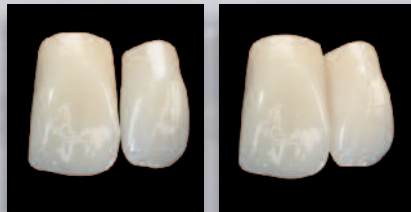

Individually layered tooth moulds

The tooth moulds are specifically layered to suit the characteristics of each individual tooth shape.



Patient Specific Moulds

In addition to the shade, the size and shape of anterior teeth should be considered in the selection process. SR PHONARES® teeth offer 18 upper and 4 lower anterior tooth moulds. Each of the moulds has been designed to deliver patient specific results.

Each of the 18 moulds are categorized by size (small, medium, large), form (soft, bold) and age (youthful, universal, mature).



SR PHONARES® *Typ* NHC

SR PHONARES®*Typ* NHC posterior teeth are suitable for a wide range of applications including partials and complete dentures.

Semi-anatomic design utilizes cusp angles of 18-22 degrees. This provides occlusion similar to anatomic teeth with decreased interference.



Conventional occlusion

Intercuspation



One-to-two setup



Quick Set Ridge lap

SR PHONARES®*Typ* teeth have been designed with a modified ridge lap that ensures minimal adjustment during setup procedures.

Full Contour

The cervical design is ideal for masking implant abutments and bars.



SR PHONARES® *Lingual* NHC

SR PHONARES® Lingual NHC posterior teeth feature progressively decreasing maxillary lingual cusps, shallow mandibular cusp angles, and an uncomplicated central fossa. They can accommodate balanced and non-balanced occlusal schemes and are suited for implant stabilized dentures.



The ideal distribution of oral forces, combined with the extraordinarily high wear resistance of the material, ensures a good balance of function and durability.


Lingualized occlusion




Intercuspatation



One-to-one setup

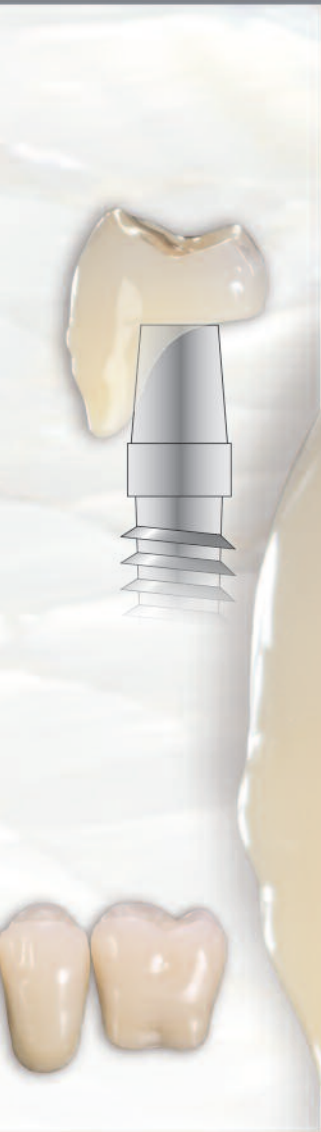


Nano Hybrid Composite for lasting results

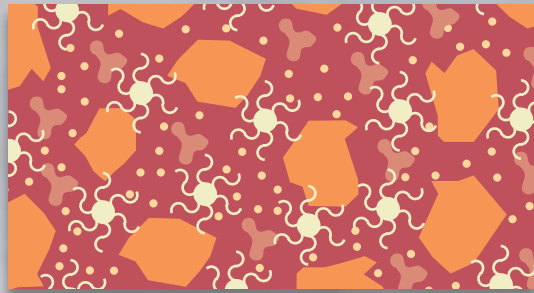


- UDMA matrix
- High-density silanized SiO₂
- Silanized SiO₂ nano-particles
- PMMA clusters polymer

In vitro wear study of denture tooth materials
Pin-on-block with steatite antagonist, after 120,000 cycles

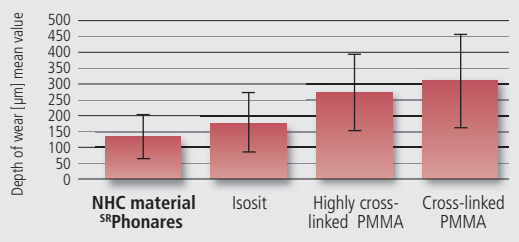

Source: Dr. Dipl. Ing. (FH) Martin Rosentritt, 08/2009, University of Regensburg, Germany

SR PHONARES® teeth are precision crafted with revolutionary Nano Hybrid Composite. This unique ISO filler technology builds on Ivoclar Vivadent's proven success with direct restoratives. The NHC material incorporates inorganic Nano fillers that deliver exceptional wear resistance, esthetics and shade stability.

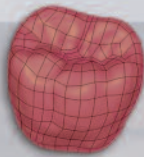

Precision CAD/CAM Design

- Precise anatomic modeling
- Improved occlusal characteristics
- Exacting posterior interdigitation




Ingenious Manufacturing

- Seamless moulding technology
- Consistent proximal contacts
- High density processing

SR PHONARES®
IMPLANT ESTHETICS

The best choice for removable implant restorations

SR PHONARES® teeth have been developed for indications when maximum wear resistance is required.

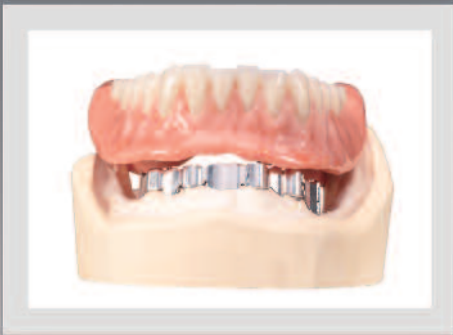

SR PHONARES® features:

- Maximum wear resistance for the most demanding indications.
- Outstanding esthetics for true to nature results
- Age specific contemporary moulds for patient specific results every time.
- Lingualized and semi anatomic occlusal schemes for both implant and conventional applications.




SHADE SELECTION

SR PHONARES® teeth are available in all 16 A-D and bleach shades for easy shade matching.



MOULD SELECTION

The anterior tooth moulds are arranged into easy to identify groups. This ensures the teeth that best match the individual characteristics of the patient can be selected.

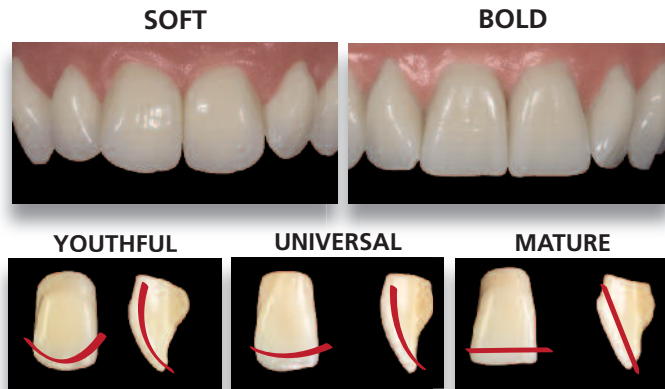

DELIVERY FORMS



SR PHONARES® NHC
18 upper anterior moulds
4 lower anterior moulds



SR PHONARES® Typ NHC
3 maxillary sets
3 mandibular sets



SR PHONARES® Lingual NHC
3 maxillary sets
3 mandibular sets

100% CUSTOMER SATISFACTION
GUARANTEED!

www.ivoclarvivadent.com

©2009 Ivoclar Vivadent, Inc.
Phonares is a registered trademark of Ivoclar Vivadent
613818 10/2009

United States

175 Pineview Drive
Amherst, NY 14228
716.691.0010
800.533.6825
Fax: 716.691.2285
mail@ivoclarvivadent.us

Canada

2785 Skymark Avenue, Unit 1
Mississauga, Ontario L4W 4Y3
905.238.5700
800.263.8182
Fax: 905.238.5711
mail@ivoclarvivadent.ca

SR PHONARES®
IMPLANT ESTHETICS

ivoclar
vivadent®
passion vision innovation